KOMPAKT-NACKEN- MASSAGEGERÄT

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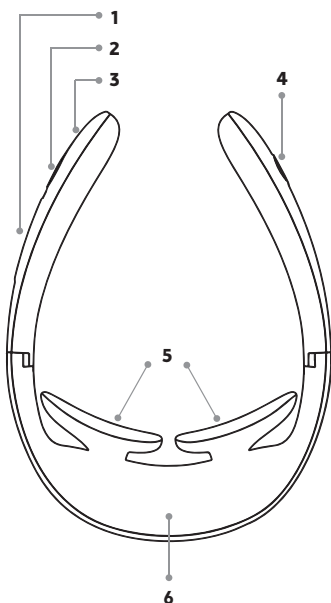
CONTENTS

At A Glance	24
Items Supplied	26
Symbols	26
Signal Words	27
Intended Use	27
When May the Device Not Be Used?	28
Safety Notices	30
Charging Up the Battery	33
Use	34
Cleaning and Storage	36
Troubleshooting	37
Disposal	39
Technical Data	39

AT A GLANCE

- 1** + / - Increase or reduce intensity
- 2**  Heat function on / off
- 3** Control lamp
 - lights up blue when the device is switched on.
 - flashes red during charging.
 - lights up red when the battery is charged.
- 4**  press for a long time: switch device on or off; press briefly: change massage type
- 5** Massage heads
- 6** USB Micro B socket (under the cover)

Dear Customer,

We normally get relaxation by warming up our muscles. Particularly for relaxing the muscles in the neck area, it may even be more pleasant to choose a mixture of muscle stimulation and heat. This provides a soothing feeling and flexibility.

We are delighted that you have chosen the **nah-vital compact neck massage device**. It can also just be placed around your neck when you are travelling and its different massage modes, heat function and 15 intensity levels provide pleasant stimulation.

If you have any questions about the device and about spare parts/accessories, contact the customer service department via our website: **www.nah-vital.de**

We hope you have a lot of fun with your **nah-vital compact neck massage device**.

Information About the Operating Instructions

Before using the device for the first time, please read through these operating instructions carefully and keep them for future reference and other users. They are an integral part of the device. The manufacturer and importer do not accept any liability if the information in these operating instructions is not complied with.

ITEMS SUPPLIED

- 1 x neck massage device
- 1 x USB cable
- 1 x operating instructions

Check the items supplied for completeness and the components for transport damage. If you find any damage, do not use the device but contact our customer service department.

Remove any possible films, stickers or transport protection from the device. **Never remove the rating plate and any possible warnings!**

A USB mains adapter is also required.

SYMBOLS

Danger symbols: These symbols indicate possible dangers. Read the associated safety notices carefully and follow them.

Supplementary information

Read the operating instructions before use!

Symbol for direct current

SIGNAL WORDS

Signal words indicate hazards if the associated notices are not followed.

DANGER – high risk, results in serious injury or death

WARNING – medium risk, may result in serious injury or death

NOTICE – may result in risk of damage to material

INTENDED USE

- This device is intended to be used for electrically stimulating the skin and the tissue underneath in the neck area.
- The device is not designed for medical or therapeutic purposes. It is not a substitute for medical treatment. If complaints persist, consult a doctor.
- The device is for personal use only and is not intended for commercial applications.
- Use the device only for the intended purpose and as described in the operating instructions. Any other use is deemed to be improper.
- The warranty does not cover faults caused by incorrect handling, damage or attempted repairs. The same applies to normal wear and tear.

WHEN MAY THE DEVICE NOT BE USED?

- Do not use the device while driving a vehicle.
- Do not use the device while you are connected to medical equipment, e.g. an ECG machine.
- Do not use the device if you are attached to high-frequency surgical devices or are in the vicinity of such devices or high-frequency commercial equipment.
- Do not use the device in the immediate vicinity of devices which operate using short waves or microwaves. Maintain a minimum gap of 1 metre.
- The device must not be used near a mobile phone. Do not use a mobile phone while using the device.
- Children, people who are sleeping or unconscious and animals must not be treated with the device.
- Do not use the device if you have a fever or physical illness/dental problems.
- Do not use the device if you have calf pain.
- The device has a heat function. People who are insensitive to heat must be careful when using the heat function.
- The following people may use the device **only following consultation with a doctor**:
 - people who are fitted with a pacemaker or similar medical implants,
 - people with heart disease or who suffer from epilepsy,
 - people who have spinal disorders or back complaints,
 - women who are pregnant or have just given birth,
 - women during their period,
 - people who suffer from circulation problems or thrombosis,

- people who have undergone an operation for which the treatment is not yet complete,
 - people with a malignant tumour,
 - people with nerve damage,
 - people under the influence of medication,
 - people who belong to a risk group or who have concerns about their health.
- Do not use the device in the vicinity of or on:
- regions which have recently been operated on or are scarred,
 - injuries and areas of swelling,
 - acne, birthmarks and warts,
 - injured, inflamed, (sun-)burnt, irritated or very sensitive skin,
 - the head, face, thyroid glands,
 - bones and joints,
 - the rib cage close to your heart,
 - the genital area,
 - body parts with no feeling.

If you are unsure whether the device may be used, consult a responsible doctor.

SAFETY NOTICES

- This device is not intended for use by people (including children) with reduced sensory or mental abilities or a lack of experience and/or knowledge, unless they are supervised by a person who is responsible for their safety or have received instruction from this person on how to use the device safely. There is a danger of injury from incorrect use!
- Do not make any modifications to the device. It may only be taken apart and / or repaired by the manufacturer, customer service department or a similarly qualified person (e.g. specialist workshop) in order to avoid any hazards.

WARNING: Do not use this device near a bathtub, shower or any vessel containing water. Never use while bathing or showering!

DANGER – Danger of Electric Shock

- Never immerse the USB cable in water or other liquids.
- If the device falls into water while it is connected to the mains power, pull the USB mains adapter out of the plug socket immediately. **Only after this** should you retrieve the device from the water.
- Never touch the device and the connecting cable with wet hands when the device is connected to the mains power.
- In the event of damage / faults, switch the device off immediately. Do not use the device if it is damaged.

WARNING – Danger from Rechargeable Batteries

- No liquid must get inside the housing as this may cause a short circuit. **Do not immerse the device in water or other liquids** and ensure that it cannot get wet!
Do not use the device if liquid has got inside it and have it checked in a specialist workshop before using it again.
- The device has a built-in battery. It must not be thrown into a fire or exposed to high temperatures of over +40°C. Do not expose it to direct sunlight for a prolonged period. Do not take apart the rechargeable battery! Do not place any heavy objects on the device. The battery could ignite.
- If the battery has leaked, do not allow your skin, eyes and mucous membranes to come into contact with the battery acid. If you do touch battery acid, rinse the affected areas with abundant clear water straight away and seek medical assistance immediately. There is a danger of burns from battery acid.

WARNING – Danger of Fire

- Keep the device away from combustible materials and substances.
- Do not cover the device during operation or charging.

WARNING – Danger of Injury

- Keep children and animals away from the packaging material. There is a danger of suffocation.
- Make sure that the USB cable is always kept out of the reach of young children and animals. There is a danger of strangulation.

- 
- The device is **not** a toy! Do not leave children and animals unsupervised with the device. There is a danger of injury from incorrect use.
 - If you experience an unpleasant feeling while using the device, stop using it immediately!
 - Stop using the device if you experience skin redness or discolouration that lasts for more than 24 hours.
 - Do not use it at the same time as other massage or stimulation devices.

Risk of Damage to Material and Property

- 
- Store the device only in closed rooms.
 - Do not expose the device to any moisture.
 - Use and store the device at room temperature.
 - Protect the device from: impacts, strong electromagnetic fields, extreme temperatures, temperature fluctuations, fire and heat sources.
 - Use only original accessories from the manufacturer in order to guarantee that there is no interference that may prevent the device from working and to avoid possible damage.
 - Disconnect the device from the power source if an error occurs during charging or before a thunderstorm.
 - Make sure that the USB cable is not squashed, bent or laid over sharp edges and does not come into contact with hot surfaces.
- 

CHARGING UP THE BATTERY

Please Note!

- Make sure that the USB cable does not present a trip hazard.
 - Do not charge up on a computer because the computer might be damaged by the high charging current.
 - Connect the USB mains adapter to a plug socket that is properly installed. It must be readily accessible so that the connection to the mains can quickly be isolated.
-

Charge up the battery fully before you first use the device, or if the device has not been used for a long time, so that it reaches full capacity.

Otherwise the battery needs to be charged if the device can no longer be switched on.

To connect the device to the plug socket, a USB mains adapter (5 V and ≥ 1 A; not supplied with the device) is required.

1. Lift up the cover of the USB socket (6).
2. First plug the USB cable into the USB socket and then connect to a power source with USB interface (preferably via a USB mains adapter to the plug socket).
The charging process begins and the control lamp flashes red. As soon as the control lamp is lit constantly red, the battery is charged.
3. Disconnect the USB cable from the power source and then from the USB socket.

USE

Beginners should use the device for no more than 15 minutes a day to prevent overstimulation on the neck!

Advanced users can use the device for up to 30 minutes a day. We recommend no more than 2 units per day, with a rest phase of at least 4 hours between the units.

Please Note!

- The device cannot be switched on while its battery is charging.
 - During the treatment, there must always be moisture between the skin and the massage heads.
 - Do not use straight after eating.
 - The device should not be used before going to sleep. Massages have a stimulating effect, which may lead to problems getting to sleep.
 - For reasons of hygiene, the device should only be used by one person and cleaned after each use (see the "Cleaning and Storage" chapter). Keep the device clean.
 - The skin may feel warm during the treatment. This is normal and not a cause for concern.
 - The massage heads must not come into contact with metal parts, e.g. on belts, wrist watches, jewellery or similar items, when the device is switched on.
 - We explicitly point out that with some people skin irritations or hypersensitivity may occur as a result of the electrical muscle stimulation. In this case, stop using the device immediately and consult a doctor if necessary.
-

1. Wipe the neck with a damp cloth.
2. Open up the clips of the device like a pair of glasses.
3. Wet the massage heads (5) with water. Use your fingers or a cotton pad to do this.
4. Hang the device around your neck so that the massage heads are resting against your neck.
5. Switch on the device. To do this, hold down the (M) button (4) for approx. 2 seconds until the device beeps once and the control lamp (3) is lit up blue.

If the device is not in contact with the skin, it switches off automatically after 3 minutes.

When the device is switched on, the intensity of the stimulation is on level 0. There are a total of 15 intensity levels.

6. Press the + button (1) repeatedly to increase the intensity to the level where stimulation still feels pleasant. Use the – button to reduce the intensity if necessary.
7. The device automatically starts the first massage mode in which four different types of massage are run through one after the other. To choose one of them permanently, briefly press the (M) button until the massage type you want starts.
After each change of massage type, the intensity is automatically set to level 0 and needs to be increased.
8. To switch the heat function on and off, press the (SSS) button (2). If it is switched on, the massage heads heat up to approx. 42°C.

The heat function can also be used separately without stimulation. To do this, press the (SSS) button immediately after switching on the device.

- 
9. The device switches off automatically after 15 minutes. To switch it off earlier, hold down the  button for approx. 2 seconds. The device beeps twice and the control lamp goes out.
 10. Clean the device after each use (see the “Cleaning and Storage” chapter).

CLEANING AND STORAGE

Please Note!

- Disconnect the device from the power supply before doing any cleaning.
 - Do not use any sharp or abrasive cleaning agents or cleaning pads to clean the device. These may damage the surfaces.
-
- Wipe down the massage heads (5) and if necessary also the other parts of the device and the USB cable with a clean damp cloth. A mild cleaning agent may also be used for more stubborn dirt.
Dry off before storing away.
 - Fold together the clips and store the device in a cool, dry place which is protected from dust and is out of the reach of children and animals.
 - If the device is not used for a prolonged period, charge up the battery every three months.

TROUBLESHOOTING

If the device does not work properly, check whether you are able to rectify the problem yourself. If the problem cannot be solved with the steps below, contact our customer service department.

Do not attempt to repair an electrical device yourself!

The device does not work.

The massage heads (5) have not been placed on correctly.

- ▶ Correct the fit. The whole of the surface of the massage heads must be in contact with the skin. There must not be any hairs, clothing, jewellery etc. in between.

There is not enough moisture between the massage heads (5) and the surface of the skin.

- ▶ Wet the skin and the massage heads with water.

The battery is too weak.

- ▶ Charge up the battery.

The stimulation is too weak.

There is not enough moisture between the massage heads (5) and the surface of the skin.

- ▶ Wet the skin and the massage heads with water.

The massage heads (5) are dirty.

- ▶ Follow the instructions in the "Cleaning and Storage" chapter.

The battery is too weak.

- ▶ Charge up the battery.

The device switches off prematurely.

The massage heads (5) have not been placed on correctly.

- ▶ Correct the fit. The whole of the surface of the massage heads must be in contact with the skin. There must not be any hairs, clothing, jewellery etc. in between.

The battery is too weak.

- ▶ Charge up the battery.

The battery is not charging. The control lamp (3) does not flash red.

The USB cable is not connected correctly to the device or the USB mains adapter.

- ▶ Check the connection and correct it.

The USB mains adapter is not inserted in the plug socket correctly.

- ▶ Correct the fit.

The plug socket is defective.

- ▶ Try another plug socket.

Pain during the treatment.

The device has been used for too long.

- ▶ Beginners should use the device for no more than 15 minutes once a day, and advanced users for no more than 15 minutes twice a day.

The skin reacts hypersensitively to the massage heads.

- ▶ Abort the treatment. If the symptoms do not subside, consult a doctor.

The massage heads are placed on a wound.

- ▶ The device must not be used on injured skin. Note the "When May the Device Not Be Used?" chapter.

DISPOSAL

Dispose of the packaging material in an environmentally friendly manner so that it can be recycled.

This device is governed by the European Directive 2012/19/EU on waste electrical and electronic equipment (WEEE). Do not dispose of the device with built-in battery as normal domestic waste, but rather in an environmentally friendly manner via an officially approved waste disposal company.

TECHNICAL DATA

Article number: 02331
 ID of operating instructions: Z 02331 M NV V2 0921 md

Device

Model number: JMK-A3
 Voltage supply: 5 V DC (1x Li-ion battery 3.7 V; 500 mAh)
 Interface: USB Micro B
 Charge time: approx. 1 h
 Operating time: approx. 30 min (with fully charged battery); automatic switch-off after 15 min
 Power: max. 3 W
 Heat function: approx. 42°C
 Mains adapter required: 5 V and ≥ 1 A

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All rights reserved.
**DS Produkte GmbH