
**FIT ENGINE**

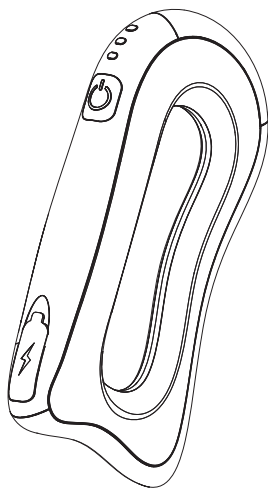**DE**

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**EN FR**

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## Contents

Items Supplied	18
At A Glance	18
Symbols	19
Signal Words	19
Intended Use	20
When May the Device Not Be Used?	20
Safety Notices	22
Charging Up the Battery	25
Use	26
Cleaning and Storage	28
Disposal	28
Troubleshooting	29
Technical Data	30

Dear Customer,

We are delighted that you have chosen the **FitEngine EMS massager**. If you have any questions about the device and about spare parts/accessories, contact the customer service department via our website: **[www.fit-engine.de](http://www.fit-engine.de)**

We hope you have a lot of fun with your **FitEngine EMS massager**.

### Information About the Operating Instructions

Before using the device for the first time, please read through these operating instructions carefully and keep them for future reference and other users. They form an integral part of the device. The manufacturer and importer do not accept any liability if the information in these operating instructions is not complied with.

EN

## Items Supplied


- EMS massage device 1 x
- USB cable 1 x
- Operating instructions 1 x

Check the items supplied for completeness and the components for transport damage. If you find any damage, do not use the device but contact our customer service department. Remove any possible films, stickers or transport protection from the device. **Never remove the rating plate and any possible warnings!**

A USB mains adapter is also required.

## At A Glance

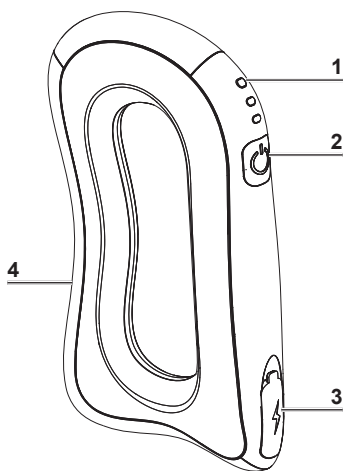- 1 3 control lamps
  - Intensity level 1, 2, 3 = blue
  - Heat function = control lamp 1 red
- 2 
  - Press for approx. 3-5 seconds: switch on / off
  - Press briefly: switch intensity and heat on / off
- 3 USB C socket (under the cover)
- 4 EMS area

## Symbols

Danger symbols: These symbols indicate possible dangers. Read the associated safety notices carefully and follow them.

Supplementary information

Read the operating instructions before use!

Direct current

## Signal Words

Signal words indicate hazards if the associated notices are not followed.

**DANGER** – high risk, results in serious injury or death

**WARNING** – medium risk, may result in serious injury or death

**CAUTION** – low risk, may result in minor or moderate injury

**NOTICE** – may result in risk of damage to material

EN

## Intended Use

- This device is intended to be used for (electrical) stimulation of muscles and fascia.
- The device is **not** designed for medical or therapeutic purposes. It is not a substitute for medical treatment. If complaints persist, consult a doctor.
- The device is for personal use only and is not intended for commercial applications.
- Use the device only for the specified purpose and as described in the operating instructions. Any other use is deemed to be improper.
- The warranty does not cover faults caused by incorrect handling, damage or attempted repairs. The same applies to normal wear and tear.

## When May the Device Not Be Used?

- Children, people who are sleeping or unconscious and animals must not be treated with the device.
- Do not use the device if you have a fever or physical illness/dental problems.
- Do not use the device while you are connected to medical equipment, e.g. an ECG machine.
- Do not use the device if you have calf pain.
- The device has a heat function. People who are insensitive to heat must be careful when using the heat function.
- The following people may use the device **only following consultation with a doctor**:
  - people who are fitted with a pacemaker or similar medical implants,
  - people with heart disease or who suffer from epilepsy,

- people who have spinal disorders or back complaints,
  - women who are pregnant or have just given birth,
  - women during their period,
  - people who suffer from circulation problems or thrombosis,
  - people who have undergone an operation for which the treatment is not yet complete,
  - people with physical abnormalities,
  - people with a malignant tumour,
  - people with nerve damage,
  - people under the influence of medication,
  - people who belong to a risk group or who have concerns about their health.
- Do not use the device in the vicinity of or on:
- regions which have recently been operated on or are scarred,
  - injuries and areas of swelling,
  - acne, birthmarks and warts,
  - injured, inflamed, (sun-)burnt, irritated or very sensitive skin,
  - bones and joints,
  - the rib cage close to your heart,
  - the genital area,
  - body parts with no feeling.

If you are unsure whether the device may be used, consult a responsible doctor.

EN

## Safety Notices

- This device is not intended to be used by people (including children) with reduced sensory or mental abilities or a lack of experience and/or knowledge, unless they are supervised by a person who is responsible for their safety or have received instruction from this person on how to use the device safely. There is a danger of injury from incorrect use!
- Do not make any modifications to the device. It may only be taken apart and / or repaired by the manufacturer, customer service department or a similarly qualified person (e.g. specialist workshop) in order to avoid any hazards.

**WARNING:** Do not use this device near a bathtub, shower or any vessel containing water. Never use while bathing or showering!

### DANGER – Danger of Electric Shock

- If the device falls into water while it is connected to the mains power, pull the USB mains adapter out of the plug socket immediately! **Only after this** should you pull the device out of the water and dispose of it.
- Never touch the device and the connecting cable (USB C cable and USB mains adapter) with wet hands when the device is connected to the mains power.

### WARNING – Danger from Rechargeable Batteries

- Make sure that the device cannot get wet. **Do not immerse in water!** Do not use the device if liquid has got inside it.

- Do not use the device if it has malfunctioned or fallen into water. Have it checked in a specialist workshop before using it again.
- The device has a built-in battery. It must not be thrown into a fire or exposed to high temperatures of over +40°C. Do not take apart the battery!
- If the battery has leaked, do not allow your skin, eyes and mucous membranes to come into contact with the battery acid. If you touch battery acid, rinse the affected areas with abundant clear water straight away and seek medical assistance immediately.

### **WARNING – Danger of Fire**

- Keep the device away from combustible materials and substances.
- Do not cover the device during operation.

### **WARNING – Danger of Injury**

- Keep the packaging material away from children and animals. There is a danger of suffocation.
- Make sure that the USB C cable is always kept out of the reach of young children and animals. There is a danger of strangulation.
- Children must be supervised to ensure that they do not play with the device. There is a danger of injury from incorrect use.
- Do not use the device while driving a vehicle.
- Do not use the device if you are attached to high-frequency surgical devices or are in the vicinity of such devices or high-frequency commercial equipment.

EN

- Do not use the device in the immediate vicinity of devices which operate using short waves or microwaves. Maintain a minimum gap of 1 metre.
- The device must not be used near a mobile phone. Do not use a mobile phone while using the device.
- If you experience an unpleasant feeling while using the device, stop using it immediately!
- Stop using the device if you experience skin redness or discolouration that lasts for more than 24 hours.
- Do not use the device if it has malfunctioned, fallen into water or if there is visible damage on the device or accessories. Have it checked in a specialist workshop before using it again.

### **Risk of Damage to Material and Property**

- Charge and store the device only in closed rooms.
- Do not expose the device to any moisture.
- Use, charge and store the device at room temperature. Do not expose it to extreme temperatures, temperature fluctuations, fire, heat sources and direct sunlight for a prolonged period.
- Protect the device from impacts.
- Do not place any heavy objects on the device or the USB C cable.
- Use only original accessories from the manufacturer in order to guarantee that there is no interference that may prevent the device from working and to avoid possible damage.
- In the event of damage / faults, switch the device off immediately.

- Do not connect the device to a plug socket during a thunderstorm.

## Charging Up the Battery

---

### Please Note!

- Do not charge up on a computer because the computer might be damaged by the high charging current.
- 

If the device cannot be switched on and control lamp 1 (1) on the device is flashing blue, the battery is too weak and needs to be charged.

To connect the device to the plug socket, a USB mains adapter (5V and  $\geq 1$ A; not supplied with the device) is required.

1. First plug the USB cable into the USB C socket (3) and then connect to a power source with USB interface (preferably via a USB mains adapter to the plug socket). The charging process begins and the control lamps on the device flash blue:
  - control lamp 1 flashes > 30%
  - control lamp 2 lights up, control lamp 2 flashes > 60%
  - control lamps 1 and 2 light up, control lamp 3 flashes > 100%
2. As soon as all three control lamps are constantly lit, the battery is charged up.  
Disconnect the USB cable from the power source and then from the USB C socket.

EN

## Use

EMS stands for electrical muscle stimulation. With this device you have a choice of treating the muscles of the neck, abdomen, legs, arms, back or buttocks.

**Beginners should treat one muscle group for no more than 15 minutes a day to prevent overstimulation!**

Advanced users can use the device for up to 30 minutes a day on each muscle group. We recommend no more than 2 units per day, with a rest phase of at least 4 hours between the units.

---

### Please Note!

- During the treatment, there must always be a little moisture between the skin and the EMS area.
- Do not use the device on your face, head or the area of the thyroid glands.
- Do not use straight after eating.
- The device should not be used before going to sleep. Massages have a stimulating effect, which may lead to problems getting to sleep.
- For reasons of hygiene, the device should be cleaned after each use (see the “Cleaning and Storage” chapter). Keep the device clean.
- The skin may feel warm during the treatment. This is normal and not a cause for concern.
- The device must not come into contact during use with metal parts, e.g. on belts, wrist watches, jewellery or similar items.
- We explicitly point out that with some people skin irritations or hypersensitivity may occur as a result of the

electrical muscle stimulation. In this case, stop using the device immediately and consult a doctor if necessary.

---

1. Clean and dry the area of the body being treated.
2. Lightly spray the area being treated with water.
3. Switch on the device:
  - Switch on / intensity level I: Press and hold  (2) for approx. 3 – 5 seconds. Control lamp 1 (1) lights up blue.
  - Intensity level II: Press  a second time (briefly), control lamps 1 and 2 light up blue
  - Intensity level III: Press  a third time (briefly), control lamps 1, 2 and 3 light up blue
  - Heat function: Press  a fourth time (briefly), control lamp 1 and the top section light up red
4. Guide the EMS area (4) over the area being treated in slow, stroking motions. Do not apply too much pressure as you do this.
5. To switch off the device, press and hold  for approx. 3 – 5 seconds. The control lamps go out.
6. Clean the device after each use (see the “Cleaning and Storage” chapter).

For more information and inspiration, have a look at our website: **[www.fit-engine.de](http://www.fit-engine.de)**.

EN

## Cleaning and Storage

### Please Note!

- Do not use any sharp or abrasive cleaning agents or cleaning pads to clean the device. These may damage the surfaces.
- Wipe down the device with a clean damp cloth. **Do not use any soap!** Dry off before storing away.
- Store in the packaging in a cool, dry place out of the reach of children and animals.

## Disposal

Dispose of the packaging material in an environmentally friendly manner so that it can be recycled.

This device is governed by the European Directive 2012/19/EU on waste electrical and electronic equipment (WEEE). Do not dispose of the device with built-in battery as normal domestic waste, but rather in an environmentally friendly manner via an officially approved waste disposal company.

## Troubleshooting

If the device does not work properly, check whether you are able to rectify the problem yourself. If the problem cannot be solved with the steps below, contact the customer service department.

**Do not attempt to repair an electrical device yourself!**

### **The device does not work.**

The EMS area (4) is not placed on correctly.

- ▶ Correct the angle or pressure slightly. There must not be any hairs, clothing, jewellery etc. between the EMS area and skin.

There is not enough moisture between the EMS area (4) and the surface of the skin.

- ▶ Moisten the skin and the EMS area with a little water.

The battery of the device is too weak.

- ▶ Charge up the battery.

### **The stimulation is too weak.**

There is not enough moisture between the EMS area (4) and the surface of the skin.

- ▶ Moisten the skin and the EMS area with a little water.

The EMS area (4) is dirty.

- ▶ Follow the instructions in the “Cleaning and Storage” chapter.

The battery of the device is too weak.

- ▶ Charge up the battery.

EN

## Technical Data

Article number: 03451  
ID of operating instructions: Z 03451 M DS V1 0721 uh

### Device

Model number: BG-S-6660  
Voltage supply: 3.7V DC (1 x Li-ion battery, 1200mAh)  
Charge time: approx. 3.5 – 4 h  
Operating time: approx. 3 h (without heat); approx. 30 – 60 min. (with heat)  
Power: 5.5W  
Heat function: 37–43°C

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All rights reserved.
DS Produkte GmbH