03775

DE

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ab Seite 3

EN

Operating instructions

starting on page 11

FR

Mode d'emploi

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NL

Bedieningshandleiding

vanaf pagina 27

Inhalt

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Contents

Symbols	11
Signal Words	11
Intended Use	12
Safety Instructions	12
Items Supplied and Device Overview	15
Functions	15
Use	16
Cleaning and Storage	17
Troubleshooting	17
Technical Data	18
Disposal	18

Symbols

Danger symbols: These symbols indicate risks of injury. Read through the associated safety instructions carefully and follow them.

Supplementary information

Read operating instructions before use!

Circuit symbol for direct current

Protection class II

Protection class III

Energy efficiency class V

May only be used indoors!

 Polarity of the coaxial connector (mains adapter)

Signal Words

Signal words indicate hazards if the associated notices are not followed.

DANGER – high risk, results in serious injury or death

WARNING – medium risk, may have serious injury or death

CAUTION – low risk, may result in minor or moderate injury

NOTICE – may result in risk of damage to material

Dear Customer,

We are delighted that you have chosen this **massage cushion**.

Before using the device for the first time, please read the operating instructions through carefully and keep them for future reference. If the device is given to someone else, these operating instructions should also be handed over with it. They form an integral part of the device.

The manufacturer and importer do not accept any liability if the information in these operating instructions is not complied with.

If you have any questions about the device and about spare parts / accessories, contact the customer service department via our website:

www.dspro.de/kundenservice

Intended Use

- This device is designed to relieve the muscles and delivers relaxation. The device can be used to massage the back, neck, abdominal, arm, foot and leg muscles. It is not a replacement for medically required treatment.
- The device is for personal use only and is not intended for commercial applications.
- The device should only be used as described in these operating instructions. Any other use is deemed to be improper.
- The following are excluded from the warranty: all defects caused by improper handling, damage or unauthorised attempts at repair. This also applies to normal wear and tear.

Safety Instructions

WARNING: Note all safety notices, instructions, illustrations and technical data provided with this device.

Failure to follow the safety notices and instructions may result in electric shock, fire and/or injuries.

- This device should not be used by people (including children) with impaired sensory or mental faculties or who do not have adequate experience and/or knowledge of such devices, unless they are supervised by a person responsible for their safety, or have received instructions from such a person on how to use the device.
- Children must be supervised to ensure that they do not play with the device.
- The device has a heating function. People who are sensitive to heat must be careful when using the heating function.

- Keep sharp and pointed objects away from the device to ensure that it is not damaged.
- The device must always be disconnected from the mains power when it is left unattended and before cleaning it.
- In the interest of your own safety, you should inspect the device for damage each time before you use it. It may only be used if it has no visible signs of damage!
- Do not make any modifications to the device. Do not replace the connection yourself. If the device or the connection are damaged, they must be replaced by the manufacturer, customer service department or a similarly qualified person (e.g. specialist workshop) in order to avoid any hazards.

Health-Specific Instructions

- **Only use the device if you have no physical illnesses or ailments.**
- The following people may use the device **only after consultation with their doctor**:
 - people who wear a pacemaker or similar medical implants;
 - people with heart disease or who suffer from epilepsy;
 - people who have back complaints or spinal disorders;
 - people who suffer from osteomyelitis, bone inflammation (in particular if the spinal column is damaged) or from osteoporosis;
 - women who are pregnant or have just given birth;
 - people who suffer from circulation problems or thrombosis;
 - people who have undergone an operation for which the treatment is not yet complete;
 - people with physical abnormalities;
 - people who have a malignant tumour.
- If you are unsure whether the device may be used, consult a responsible doctor.
- Do not use the device near to or on injuries, swelling, skin irritations, infected skin, warts, fungal infections, abnormal changes in the skin or body parts with no feeling.
- The device must not be used for longer than 30 minutes continuously. Do not massage one group of muscles for longer than 15 minutes.
- If you experience an unpleasant sensation while you are using the device, stop performing the massage immediately!
- Do **not** fall asleep while you are using the device!
- Do not use the device before going to bed. The massage has a stimulating effect, which may lead to problems getting to sleep.

**DANGER – Danger of Electric Shock**

- Use and store the device only in closed rooms.
- Do not operate or store the device in rooms with high humidity.
- Never immerse the device or the connection in water or other liquids and ensure that they cannot fall into water or become wet.
- Should the device fall into water, switch off the power supply immediately. Do not attempt to pull the device out of the water while it is still connected to the mains power!
- Never touch the device or the connection with wet hands when the device is connected to the mains power.

**WARNING – Danger of Fire**

- Do not use the device under a blanket or similar item or in the vicinity of combustible material.

**WARNING – Danger of Injury**

- Keep plastic film and bags well away from children and animals. There is a danger of suffocation.
- Make sure that the connection is always kept out of the reach of small children and animals. There is a danger of strangulation.
- Do not insert any fingers or other body parts into the device. Do not open the zip during use.
- Make sure that the connected mains cable does not present a trip hazard. There is a danger of falling.


NOTICE – Risk of Damage to Material and Property

- Only connect the mains adapter to a plug socket that is properly installed and matches the technical data of the device. The plug socket must also be readily accessible after connection so that the connection to the mains can quickly be isolated.
- Only use suitable extension cables whose technical data is the same as that of the device.
- In the event of damage / faults (e.g. unusual noises, odour or smoke development) and in a thunderstorm, switch the device off immediately and pull the mains adapter out of the plug socket.
- The device must be switched off when inserting or removing the mains adapter. Always pull the mains adapter and never the mains cable to disconnect the device from the mains power.
- Do not lie or stand on the device.
- Protect the device from heat, naked flames, direct sunlight, extreme temperatures, moisture and collisions.
- Keep the device in a dry place at temperatures of between 0 and 40°C.

Items Supplied and Device Overview

Check the scope of delivery for completeness and the components for transport damage. If you find any damage, do not use the product but contact our customer service department.

Remove any possible films, stickers or transport protection from the product. **Never remove the rating plate and any possible warnings!**

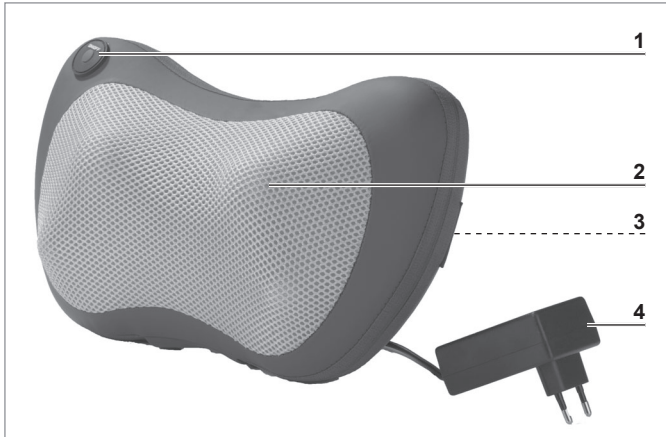- 1 **ON/OFF** button
- 2 2 x 2 massage balls, heating and red light element
- 3 Fastening elastic band
- 4 Mains adapter

Functions

You can use the **ON/OFF** button (1) to choose between the following options:

- Massage function with automatic change of direction every 5 minutes, red light and heat, starts with clockwise rotation (press the **ON/OFF** button once)
- Massage function with automatic change of direction every 5 minutes, red light and heat, starts with anticlockwise rotation (press the **ON/OFF** button twice)
- Massage function with automatic change of direction every 5 minutes, without red light and heat (press the **ON/OFF** button three times)
- Switch off (press the **ON/OFF** button four times)

Use

PLEASE NOTE!

- No pressure should be exerted on the device as this will slow down or stop the massage and the device could be damaged. Do not stand or lie on the device.
- Do not use the device for longer than 30 minutes without stopping.
- Do not massage one group of muscles for longer than 15 minutes.
- For reasons of hygiene, avoid direct contact of the skin with the device wherever possible. It should not come into contact with skin that has had cream or oil applied to it.
- A tingling sensation/itching during or after the massage is normal due to the circulation-boosting effect.

We recommend that you wait an hour after mealtimes before using the device.

Areas of Application

When this device performs a massage, two massage heads massage the muscles twice. This is intended to relax them and boost circulation. The device is suitable for massaging the back, neck, abdominal, arm, foot and leg muscles.

Examples:

Neck and shoulders

Back

Thighs

Operating the Device

1. Connect the mains adapter to the connecting cable on the device.
2. Insert the mains adapter into an easily accessible socket.
3. Place the massage balls (2) of the device onto the part of the body which is to be massaged (see the "Areas of Application" section). For a shoulder or back massage, pull the fastening elastic band (3) over the back of a chair. Make sure it is the correct height.
4. Press the **ON/OFF** button (1) several times if necessary (see the "Functions" chapter). The massage starts.
5. After 15 minutes the device automatically switches off. To switch it off sooner, press the **ON/OFF** button (1) until the device switches off.
6. Pull the mains adapter out of the plug socket and clean the device if necessary.

Cleaning and Storage

PLEASE NOTE!

- The device may not be connected to the mains and must be completely switched off for cleaning!
 - Do not use any caustic or abrasive cleaning agents for cleaning. These can damage the surface.
-
- Wipe down the device just with a damp cloth if necessary. Then wipe it down fully with a dry cloth.
 - Do not wind the mains cable around the device or the mains adapter.
 - Keep the device in its respective packaging in a dry place which is protected from sunlight and is not accessible to children and animals.

Troubleshooting

If the device does not work properly, check whether you are able to rectify the problem yourself. If the problem cannot be solved with the steps below, contact the customer service department.

Do not attempt to repair a electrical device yourself!

Problem	Possible cause / solution
The device does not work.	<ul style="list-style-type: none"> • Has the ON/OFF button (1) been pressed properly? • Is the mains adapter inserted in the plug socket? • Is the plug socket defective? Try another one. • Check the fuse of the mains connection. • Is the connecting cable connected correctly to the mains cable of the mains adapter? • The device has switched off automatically after operating for 15 minutes. Press the ON / OFF button (1) again. • There is too much weight on the device.
No heat is produced.	<ul style="list-style-type: none"> • Is the red light and heat switched on?
The device gets too warm.	<ul style="list-style-type: none"> • Was it in operation several times in succession for 15 minutes / too long in each case? Pull the mains adapter out of the plug socket and allow the device to cool down. Always disconnect the device from the mains when it is not in use.
The device makes noises.	<ul style="list-style-type: none"> • Operating noises are a normal consequence of the mechanics and not a defect.
Smoke emerges from the device.	<ul style="list-style-type: none"> • Pull the mains adapter out of the plug socket immediately and contact the customer service department.

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All rights reserved.

Disposal

Dispose of the packaging material in an environmentally friendly manner so that it can be recycled.

This device is governed by the European Directive 2012/19/EU on waste electrical and electronic equipment (WEEE). Do not dispose of this device as normal domestic waste, but rather in an environmentally friendly manner via an officially approved waste disposal company.

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Gefabriceerd voor DS Produkte GmbH