INHALT

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CONTENTS

Items Supplied	10
At A Glance	10
Symbols	11
Signal Words	11
Intended Use	12
When May You Not Use The Device?	12
Safety Notices	13
Before Initial Use	15
Battery	15
Attachments and Areas of Application	16
Use	16
Care and Storage	16
Troubleshooting	17
Disposal	17
Technical Data	18

ITEMS SUPPLIED

(Picture A)

- A Fork attachment
- B Flat attachment
- C Cartridge attachment
- D Sphere attachment
- E Main unit
- F Mains adapter
- Storage case
- Operating instructions

Check the items supplied for completeness and the components for transport damage. If you find any damage, do not use the device but contact our customer service department.

AT A GLANCE

(Picture B)

- 1 Intensity level control lamps
- 2  Start / stop / set intensity level
- 3 Operating lamp
- 4 On/Off switch
- 5 Connecting socket for mains adapter
- 6 Battery control lamp
- 7 Handle
- 8 Mount for the attachment
- 9 Ventilation slots

Dear Customer,

We are delighted that you have chosen our **MAXXMEE intensive massage device**. This device can help you to shorten the recovery time after peak levels of performance, promote blood flow and restore your flexibility.

If you have any questions about the device and about spare parts/accessories, contact the customer service department via our website:

www.dspro.de/kundenservice

We hope you have a lot of fun with your **MAXXMEE intensive massage device**.

Information About the Operating Instructions

Before using the device for the first time, please read through these operating instructions carefully and keep them for future reference and other users. They form an integral part of the device. The manufacturer and importer do not accept any liability if the information in these operating instructions is not complied with.

SYMBOLS

Danger symbols: These symbols indicate possible dangers.

Read the associated safety notices carefully and follow them.

Supplementary information


Read operating instructions before use!

Symbol for direct current

Only use indoors!

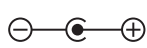
Polarity of the coaxial connector (mains adapter)

Protection class II

Safe electrical isolation

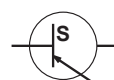Switching power supply

Energy efficiency class VI

SIGNAL WORDS

DANGER

warns of serious injuries and danger to life

WARNING

warns of **possible** serious injuries and danger to life

NOTICE

warns of material damage

INTENDED USE

- This device is intended to be used for massaging the whole body and delivers relaxation. **It is not a replacement for medically required treatment.**
- Do not use the device before going to bed. Massages have a stimulating effect, which may lead to problems getting to sleep.
- The device is for personal use only and is not intended for commercial applications.
- The device should only be used as described in the operating instructions. Any other use is deemed to be improper.
- The warranty does not cover faults caused by incorrect handling, damage or attempted repairs. The same applies to normal wear and tear.

WHEN MAY YOU NOT USE THE DEVICE?

- The following people may use the device **only after consultation with their doctor**:
 - people who wear a pacemaker or similar medical implants;
 - people with heart disease or who suffer from epilepsy;
 - people who have back complaints or spinal disorders;
 - people who suffer from osteomyelitis, bone inflammation (in particular if the spinal column is damaged) or from osteoporosis;
 - women who are pregnant or have just given birth;
 - people who suffer from circulation problems or thrombosis;
 - people who have undergone an operation for which the treatment is not yet complete;
 - people with physical abnormalities;
 - people who have a malignant tumour.
- Do not use the device near to or on:
 - injuries,
 - any swelling,
 - skin irritation,
 - infected skin,
 - warts,
 - fungal infections,
 - abnormal changes in the skin,
 - bones,
 - head,
 - genital area,
 - body parts which are numb or have no feeling.
- **Only people who have been trained may use the device on joints or all around joints.**

If you are unsure whether the device may be used, consult a responsible doctor.

SAFETY NOTICES

WARNING: Note all safety notices, instructions, illustrations and technical data provided with this device. Failure to follow the safety notices and instructions may result in electric shock, fire and/or injuries.

- This device should not be used by people (including children) with impaired sensory or mental faculties or who do not have adequate experience and/or knowledge of such devices, unless they are supervised by a person responsible for their safety, or have received instructions from such a person on how to use the device.
- **Children** must not play with the device.
- Cleaning and **user maintenance** may not be performed by **children**, unless they are supervised.
- **Children** younger than 8 years of age and animals must be kept away from the device and connecting cable.
- Sleeping or unconscious people and animals may not be treated with the device.
- The device must always be disconnected from the mains power when left unattended and before cleaning it.
- Do not make any modifications to the device. If the device or the connecting cable are damaged, they must be repaired/replaced by the manufacturer, customer service department or a specialist workshop in order to avoid any hazards.

DANGER – Danger of Electric Shock

- Use, charge and store the device only in closed rooms.
- Do not use the device in rooms with high humidity.
- Never immerse the device in water or other liquids! Ensure that the device and the connecting cable cannot fall into water or become wet.
- If the device does fall into water, pull the mains plug out of the plug socket immediately. Do not attempt to pull the device out of the water while it is still connected to the mains power!
- Never touch the device and the connecting cable with wet hands when the device is connected to the mains power.
- Do not use the device if it has malfunctioned or been dropped in water. Have it checked in a specialist workshop before using it again.

WARNING – Danger of Injury

- **Danger of suffocation!** Keep children and animals away from the packaging material.
- **Danger of strangulation!** Make sure that the connecting cable is always kept out of the reach of small children and animals.
- **Danger of burns from battery acid!** If the battery leaks, avoid physical contact with the battery acid. Wear protective gloves if necessary. If you touch battery acid, rinse the affected areas with abundant clear water straight away and seek medical assistance immediately.
- In the interest of your own safety, inspect the device and the connecting cable for damage each time before you use it. Use the device only if it is in working order.

- **Danger of tripping!** Make sure that the connecting cable does not present a trip hazard.

WARNING – Danger from Rechargeable Batteries

- The main unit with a built-in battery must not be taken apart, thrown into a fire or exposed to high temperatures of over +40°C. The battery could ignite or explode.
- Make sure that the device, and in particular the main unit, cannot get wet. Do not use the device if liquid has got into it. The battery in the main unit could ignite or explode.
- Do not cover the device and the mains adapter during charging in order to prevent any overheating.

NOTICE – Risk of Damage to Material and Property

- The device must be switched off when the mains adapter is pulled out or inserted into the plug socket.
- Make sure that the connecting cable is not squashed, bent or laid over sharp edges and does not come into contact with hot surfaces.
- When the mains adapter is removed from the plug socket, always pull on the mains adapter and never the mains cable.
- Do not pull or carry the device by the mains cable. Never wrap the mains cable around the mains adapter.
- Do not insert anything into the ventilation openings of the device and make sure that these do not become clogged. The device could otherwise overheat and be damaged.
- Use only original accessories from the manufacturer in order to guarantee that there is no interference that may prevent the device from working and to avoid possible damage.

- Protect the device from heat sources, naked flames, sub-zero temperatures, persistent moisture, wet conditions and impacts.

BEFORE INITIAL USE

- Remove any possible films, stickers or transport protection from the device. **Never remove the rating plate and any possible warnings!**
- Charge up the battery (see the "Battery" chapter – "Charging Up the Battery").

BATTERY

Charging Up the Battery

Please Note!

- **Risk of short-circuits!** Do not insert any conducting objects (e. g. screw-drivers or similar items) into the charging socket.
- Only connect the mains adapter to a plug socket that is properly installed and matches its technical data. The plug socket must be readily accessible after connection so that the connection to the mains can quickly be isolated.
- Only use suitable extension cables whose technical data is the same as that of the mains adapter.
- Unwind the connecting cable fully before connecting it.
- Remove the mains adapter from the plug socket if a fault occurs during charging or if a thunderstorm is expected.
- Charge the battery only with the mains adapter supplied with it. Do not use the mains adapter that is supplied to charge up any other devices.

- **It is essential that you charge up the battery fully before you first use the device in order to achieve the maximum charging capacity.**

If the battery control lamp (6) and the operating lamp (3) are lit red, the battery is weak and should be charged.

1. Make sure that the device has been switched off. To do this, move the On/Off switch (4) to the left.
2. First plug the barrel connector of the mains adapter (F) into the connecting socket (5) of the main unit (E) and then connect the mains adapter to a plug socket. The battery control lamp flashes green during the charging process.
The battery is fully charged when the battery control lamp is lit up constantly green.
3. First disconnect the mains adapter from the mains power and then disconnect the barrel connector from the device.

Battery Status Displays

The battery control lamp (6) and the operating lamp (3) indicate the battery status by their colour:

- **green:** battery full
- **amber:** battery half full
- **red:** battery almost empty

NOTICE! The battery must not be fully discharged as this can cause damage to the battery! To protect the battery, the device therefore switches off automatically when the battery is weak after the battery control lamp and the operating lamp have flashed twice. If this happens, the battery should be charged!

ATTACHMENTS AND AREAS OF APPLICATION

- **Sphere attachment** (Pictures C): for large muscle groups such as quadriceps and gluteus maximus and for targeted loosening if you have aching muscles
- **Flat attachment** (Pictures D): suitable for relaxing the whole body
- **Cartridge attachment** (Pictures E): for deeper tissue and for the feet, hands etc.
- **Fork attachment** (Pictures F): for the neck, back extensor, calves or similar muscles.

USE

Please Note!

- To prevent overheating, do not use the device under a blanket or similar cover.
- Do not use the device while it is connected to the mains power.
- For reasons of hygiene, the sphere attachment should not be used on skin that has had cream or oil applied to it.
- If you experience a feeling of discomfort or tingling in the hand you are holding the device with, stop using the device and consult a doctor if necessary.

Check the battery status before each use (see the "Battery" chapter – "Battery Status Displays").

1. Insert the attachment you want as far as it will go into the mount (8) on the main unit (E).
2. Move the On/Off switch (4) to the right to switch on the device. The battery control lamp (6) and the operating lamp (3) light up. The device is in standby mode.

3. Press the button  (2) to start the device and set the intensity level:
 - Press 1x: Intensity level 1
 - Press 2x: Intensity level 2
 - Press 3x: Intensity level 3
 - Press 4x: Off
 The number of lit control lamps (1) indicates which intensity level has currently been set.

We recommend starting with the lowest level and gradually increasing the intensity.

4. Hold the device by the handle (7) and guide the attachment over the area of the body you want, applying slight pressure. **Do not apply too much pressure!**

We recommend that you do not massage one muscle area for longer than a maximum of 15 minutes at a time.

CARE AND STORAGE

Please Note!

- Do not use any sharp or abrasive cleaning agents or cleaning pads to clean the device. These may damage the surfaces.

Make sure that the device is switched off and is not connected to the mains power.

- Remove the used attachment from the main unit (E) and clean it. All attachments **apart from the sphere attachment** (D) can be cleaned with warm water and a mild detergent.
- Wipe down the main unit and the sphere attachment with a slightly damp cloth.

- Store the device and all accessories in a storage case in a dry place that is protected from sunlight and is not accessible to children and animals.

TROUBLESHOOTING

If the device does not work properly, check whether you are able to rectify the problem yourself. If the problem cannot be solved with the steps below, contact the customer service department.

Do not attempt to repair an electrical device yourself!

The battery control lamp (6) does not flash when the device is connected to the mains power.

The mains adapter (F) is not inserted correctly in the plug socket or the barrel connector is not inserted correctly in the connecting socket (5) of the device.

- ▶ Correct the fit of the plug-in connections.

The plug socket is defective.

- ▶ Try another plug socket.

There is no mains voltage.

- ▶ Check the fuse of the mains connection.

The device does not start.

The battery is empty.

- ▶ Charge up the battery.

The On/Off switch (4) has been moved to the left.

- ▶ Move the On/Off switch to the right.

The button  (2) has not been pressed.

- ▶ Press the  button to switch on the device.

DISPOSAL

Dispose of the packaging material in an environmentally friendly manner so that it can be recycled.

This device is governed by the European Directive 2012/19/EU on waste electrical and electronic equipment (WEEE). Do not dispose of the device with built-in battery as normal domestic waste, but rather in an environmentally friendly manner via an officially approved waste disposal company.

TECHNICAL DATA

Article number: 07770
ID of operating instructions: Z 07770 M DS V1 0320 dk

Device

Model number: 72841A0
Power supply: 22.2 V (Li-ion battery, 2400mAh)
Intensity levels: 3 levels (1 = 1800 strokes / min.; 2 = 2400 strokes / min.; 3 = 3200 strokes / min.)
Battery charge time: approx. 5 hours
Battery operating time: up to approx. 3 hours (when the battery is fully charged)

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Modelnummer: GQ24-260100-AG