SPORTS ONLY X VENENWALKER®

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EN Operating instructions

ES Manual de instrucciones

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DE ab Seite 3

EN from page 20

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NL vanaf pagina 54

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IT da pagina 88

EMC Manufacturer's Declaration

Seite / page / pagina / página 105

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XIAMEN SENYANG CO.,LTD

Adresse: No 189-1, Tongji South Road, Jimei District, Xiamen, 361000, Fujian, P.R. China

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LIST OF CONTENTS

Control unit at a glance **21**

Symbols on the Device and in the Operating Instructions **22**

Signal Words **22**

Purpose **24**

Brief description **24**

Use **24**

Safety Notices **26**

Items Supplied **28**

Before Initial Use **28**

Use **29**

Cleaning **32**

Storage **32**

Troubleshooting **33**

Technical Data **34**

Customer Service **35**

More Information **35**

Disposal **36**

IMPORTANT NOTICES! **PLEASE KEEP FOR REFERENCE!**

Before using the device for the first time, please read through these operating instructions carefully and keep them for future reference and for other users. They form an integral part of the device.

The manufacturer and importer do not accept any liability if the information in these operating instructions is not complied with.


The illustrations in these operating instructions may differ from your device.

CONTROL UNIT AT A GLANCE

- 1 Handle (can be pulled out)
- 2 Socket for connecting the connection tubes
- 3 Control panel

Control panel

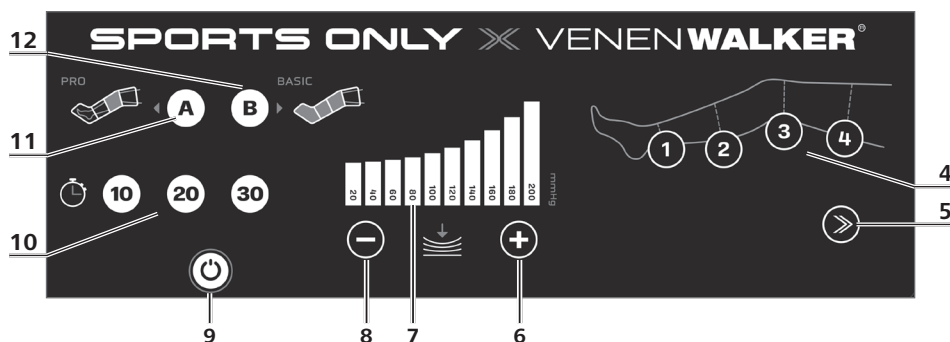- 4 **① – ④** Selection and display of chamber activity: When you touch one of the buttons, the associated chamber is activated (button is lit) and deactivated (button is not lit).
- 5 **➤** Start/stop: blue = device is stopped; green = device is operating
- 6 **+** Increase cuff pressure
- 7 Display of the chosen cuff pressure
- 8 **–** Reduce cuff pressure
- 9 **⏻** Main power switch: Touch briefly to switch on the device; hold for a few seconds to switch it off.
- 10 **10 20 30** Timer: switch off automatically after 10, 20 or 30 minutes
- 11 **A B** Mode **A** (Pro mode): Pressure switches from one chamber to the next (recommendation: sports application).
- 12 **A B** Mode **B** (Basic mode): Pressure is built up from the foot upwards and is maintained in the lower chambers (recommendation: medical application).

SYMBOLS ON THE DEVICE AND IN THE OPERATING INSTRUCTIONS

Danger symbols: These symbols indicate hazards. Read through the associated safety notices carefully and follow them.

Supplementary information

Read the operating instructions before use! This symbol appears in blue on the device.

Do not wash

Do not bleach

Do not dry in a tumble dryer

Do not iron

Do not dry-clean

Protection class II

Symbol for devices classified as "Type BF" (level of protection against electric shock)

Manufacturer

Date of manufacture

Serial number (date of manufacture / sequential number)

Batch code

Alternating current

Devices bearing this symbol were put into circulation after 13.08.2005 (DIN EN 50419). Devices labelled with this symbol meet the labelling requirements of all EU member states.

International recycling symbol

Temperature range

Humidity range

Air pressure range

Keep dry

Fragile

EU Representative

Importer

Medical product

Sales

WITH THE CE MARK, THE MANUFACTURER OR IMPORTER CONFIRMS THAT THE PRODUCT MEETS THE REQUIREMENTS OF THE APPLICABLE EC DIRECTIVES AND COMPLIES WITH THE "GENERAL REQUIREMENTS" STIPULATED THEREIN. The number 0197 alludes to the notified body (certifier). Together these certify that the Medical Product Law Implementing Act (MPDG, Medizinproduktegesetz) has been complied with by the distributor.

SIGNAL WORDS

Signal words indicate hazards if the associated notices are not followed.

DANGER – high risk, results in serious injury or death

WARNING – medium risk, may result in serious injury or death

NOTICE – may result in risk of damage to material

Dear Customer,

We are delighted that you have chosen **Sports Only from VenenWalker®**.

This device will help your legs to recover after intensive training. If you have a medical indication, please seek medical advice before using the device.

Sports Only from VenenWalker® is a medical product certified according to DIN EN ISO 13485:2016 and meets the necessary safety and performance requirements according to MDR:2017: (Appendix I).

If you have any questions about the device and about spare parts/accessories, contact the customer service department via our website: **www.venenwalker.com** or by phone on +49 (0) 38851 314337* (Mon – Fri, 8 am – midday)*

We wish you every success using the **Sports Only from VenenWalker®!**

Looking After Yourself

VenenWalker® has been a well-known brand for medical devices that operate using intermittent (alternating) compression since 1997. The use of such devices in Germany is described in detail in the guidelines for intermittent compression devices. According to these guidelines, the cuff pressure must not exceed 120 mmHg.

This restriction does not apply in some other countries. These countries also permit the use of medical devices that produce a pressure of 240 mmHg, and in some cases much higher. In our assessment, this is because these countries do not have very much experience in the field of intermittent compression.

As our certified medical product **Sports Only from VenenWalker®** can also be used to produce a pressure higher than 120 mmHg, the device should only be used, if a **pressure of over 120 mmHg is set, by people with healthy legs and no medical indication.**

Use the device **only at a pressure of no more than 120 mmHg if you have any medical indications!** A higher pressure could cause health problems. However, in the range up to and including 120 mmHg, the device can be used safely even if you have medical indications.

IMPORTANT NOTICE FOR THE USERS OF THIS MEDICAL PRODUCT

Please notify the manufacturer, the EU Representative or the importer about any serious incidents that have occurred in relation to the medical product!

A serious incident is an incident that causes serious limitation to people's health.

In particular, serious incidents include any incidents that may have resulted directly or indirectly in death or do result in death or cause a temporary or permanent serious deterioration in the health condition of a patient, user or another person or represent a serious threat to public health.

* Calls to German landlines are subject to your provider's charges.

PURPOSE

- **Sports Only from VenenWalker®** (also referred to hereinafter as the "device") may only be used for decongestion of the extremities. **If you have a medical indication, please seek medical advice before using the device!**
- The device is for personal use only and is not intended for commercial applications.
- Use the device only for the intended purpose and as described in the operating instructions. Any other use is deemed to be improper.
- The warranty does not cover faults caused by incorrect handling, damage or attempted repairs. The same applies to normal wear and tear.

BRIEF DESCRIPTION

Sports Only from VenenWalker® is a system comprising a control unit and two double-walled leg cuffs with four interior chambers which also cover the feet and soles.

The control unit feeds air into the chambers of the cuffs via connection tubes that are connected to it. This enables controlled pressure to be applied so as to pump up the cuffs gently and/or massage the legs.

In addition to its function, with the **Sports Only from VenenWalker®**, as with all our products, we attach particular importance to user safety. **Please note that a pressure of over 120 mmHg may have an adverse effect on the arterial blood supply.**

USE

The device is designed in such a way that no special knowledge or skills are required to use it. It may be used up to three times per day for 10 to 30 minutes on each occasion. The device is suitable for long-term use.

The indications and contraindications listed below are taken from the current guideline for apparative intermittent compression (AIC or IPC) issued by the German Society of Phlebology. For more detailed information, please consult the "S1 Guideline 037-007, Intermittent Pneumatic Compression (IPC, AIC)" version 20.06.2017 incl. editorial revision SR 11.01.2018, published by AWMF online – the scientific medical platform.

Indications

The apparative intermittent compression (AIC) of the **Sports Only from VenenWalker®** is intended to be used exclusively:

- for recovery after intensive training
- for warm-up (activation) before doing sport
- to treat congestion caused by joint injury or overexertion, e.g. haematoma and sore muscles in the event of sports injuries
- for medical indications up to a pressure of 120 mmHg

Contraindications

If you suffer from the following complaints, we recommend the **VenenWalker Pro2**, the **VenenWalker basic**, or the **Sports Only up to 120 mmHg**, but **NOT** the **Sports Only from VenenWalker®** with a pressure of over 120 mmHg:

- Leg swelling caused by lack of movement, for example due to standing for long periods of time in the case of sales personnel, due to sitting for extended periods when doing office work or due to long-haul travel
- Swollen legs due to being bedridden or as a result of pregnancy**
- Swollen legs following a medical diagnosis
- Venous diseases such as varicose veins and spider veins
- Lipoedemas
- Prevention of embolisms due to being bedridden
- Lymphatic system swelling (lymphoedema)
- Lower-leg ulcers caused by a venous circulation disorder
- Congestion due to joint injury or overexertion, e.g. haematoma and sore muscles in the event of sport injuries
- Swelling following a pulled muscle or sprain
- To assist the wound healing process after an operation
- Peripheral arterial occlusive disease, with certain indications**

Contraindications for Sports Only from VenenWalker®, VenenWalker® Pro2 and VenenWalker® basic

In the following cases, do **NOT** use the **Sports Only from VenenWalker®, VenenWalker® Pro2 or VenenWalker® basic**, or do so only after consulting a doctor:

- Decompensated cardiac insufficiency
- Extended thrombophlebitis, thrombosis or suspicion of thrombosis
- Acute erysipelas
- Severe, unmanaged hypertension
- Pronounced neuropathy of the extremities
- Acute phlegmons
- Compartment syndrome
- Extended or open soft-tissue trauma of the extremities
- Occluding processes in the lymph drainage area where AIC has caused congestion in the groin and genital area.
- Blistering dermatoses
- Pregnant women may only use the device after consulting a doctor.

If you are not sure whether you suffer from one of the above illnesses, seek medical advice before using the system!

** Seek medical advice.

Precautionary measures

The **Sports Only from VenenWalker®** should be used with care on patients with the following symptoms or conditions:

- Peripheral neuropathy, pain or numbness of the limbs
- Undiagnosed, untreated or infected wounds, inflamed skin, transplantations or dermatological conditions which are exacerbated by the cuff
- Extremely deformed limbs which prevent the correct application of the cuff

SAFETY NOTICES

- This device should not be used by people (including **children**) with impaired sensory or mental faculties or who do not have adequate experience and/or knowledge of such devices, unless they are supervised by a person responsible for their safety, or have received instructions from such a person on how to use the device.
- Children must **not** play with the device.
- The cleaning may not be carried out by **children**.
- Children under 8 years of age and animals must be supervised near the device.
- Electrical and medical devices can be dangerous if used improperly. The control unit does not have any parts which can be serviced by the user. The housing of the control unit must not be opened. Under no circumstances may the control unit be modified!
- If components of the device are damaged, they must be replaced or repaired by the manufacturer, customer service department or authorised technical personnel in order to avoid hazards.
- The control unit may only be used with accessories approved by the manufacturer, and the accessories may not be used with other compressors. Otherwise it is not possible to guarantee that the device will work correctly, and all guarantee and warranty claims will be invalid.
- The device must always be disconnected from the mains power when left unattended and before cleaning it.

DANGER – Danger of Electric Shock

- Use and store the device only in closed rooms.
- Do not operate or store the device in rooms with high humidity.
- Never immerse the control unit, its connecting cable and the connection tubes in water or other liquids and ensure that these components cannot fall into water or become wet.
- If the control unit falls into water while it is connected to the mains power, switch off the power supply immediately. **Only after this** should you pull the device out of the water! Have the control unit checked in a specialist workshop before using it again.
- Never touch the control unit or its connecting cable with wet hands when the control unit is connected to the mains power.

WARNING – Danger of Injury

- Keep the packaging material away from children and animals. There is a danger of suffocation.
- Make sure that the connecting cable and the connection tubes are always kept out of the reach of small children and animals. There is a danger of strangulation.
- Make sure that the connecting cable and the connection tubes do not present a trip hazard.
- Sleeping or unconscious people may **not** be treated with the device.
- Treatment should be interrupted immediately if pain, a tingling sensation or numbness occurs in the limbs during or after treatment.
- In the event of a fault in which the cuff remains inflated, disconnect the connection tubes so as to let the air out of the cuff(s) and remove the cuff(s) from the limbs.
- Do not use the control unit if it has malfunctioned. Have it checked in a specialist workshop before using it again.

NOTICE – Risk of Damage to Material and Property

- Keep the cuffs and connection tubes away from sharp objects. Do not insert safety pins or other pins or sharp objects into the cuffs.
- Make sure that the mains cable, the cuffs and the connection tubes are not crushed, bent or laid over sharp edges and do not come into contact with hot surfaces.
- Switch off the device and pull the mains plug out of the plug socket if an error occurs during operation or before a thunderstorm.
- When removing the mains plug from the plug socket, always pull on the mains plug and never the mains cable. Do not pull or carry the control unit by the connecting cable.
- Keep the control unit away from highly flammable materials and do not cover it during operation to prevent overheating or a fire. Do not insert anything into the ventilation openings of the control unit and make sure that these do not become clogged.
- Always place the control unit on a dry, even and firm surface. Do not place on cushions or other soft surfaces in order to prevent the control unit from overheating.

- Do not place the device on or next to heat sources such as a radiator, hobs or ovens. Keep flames away from the device.
- Protect the device from flames, extreme temperatures, persistent direct sunlight and impacts. Please refer to the chapter entitled "Technical Data" for details.
- Never place heavy objects on the control unit, the cuffs or connection tubes.

Electromagnetic Compatibility (EMC)

For details, please refer to the accompanying documents at the end of these operating instructions (EMC Manufacturer's Declaration) or go to

<https://www.venenwalker.com/downloads/>

ITEMS SUPPLIED

- 1 x control unit
- 2 x cuffs, each with 4 connection tubes
- 1 x rubber foot
- 1 x operating instructions

- Carefully take all parts out of the packaging and remove any films, but **never remove the rating plate and any possible warnings!** Do not use any scissors or other sharp or pointed objects to unpack the cuffs so you do not damage them.


We recommend that you keep the packaging so you can store the device in it safely when it is not in use.

- Check the items supplied for completeness and the components for transport damage. If you notice any damage, do not use the device, but contact the customer service department.

BEFORE INITIAL USE

After unpacking the control unit, remove the transport safeguard as follows:

- On the bottom of the control unit there is a wing screw in the middle. Unscrew the wing screw (1.) and replace it with the rubber foot that is also supplied (2.).

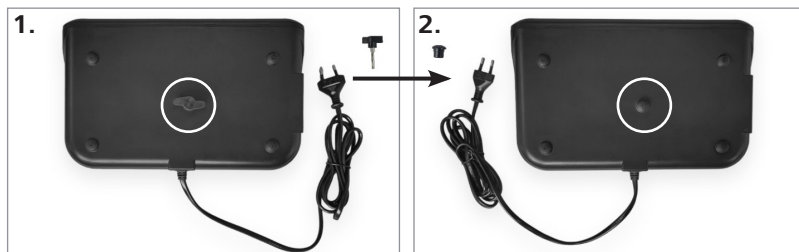USE

Mode of Operation

The device works based on apparative intermittent compression (AIC). Your legs are massaged by the rhythmic, alternating inflation and deflation of the air cushions. The cuffs exert a regularly interrupted pressure from outside on the capillary pump in the calf muscles and feet. The water inside the tissue and the metabolic particles dissolved in it are flushed out and can be eliminated from the body.

Notices

- Your legs should be horizontal or raised when using the system so as to allow the blood to circulate unhindered. Assume a horizontal position if possible.
- We recommend using the device for approx. 10 to 30 minutes in the morning and evening.
- Do not use the device for long periods without interruption since excessive massage could be too strenuous. When you use the device for the first time, limit yourself to 10 minutes. If you are comfortable with a period of 10 minutes, you can extend the length of subsequent treatments.

A) Setting up the Control Unit

PLEASE NOTE!

- Always check the device for any damage or wear before putting it into operation. It may only be used if it is undamaged and does not display any signs of wear and tear.
 - If the control unit is moved from a very cold environment to a warm environment, condensation may collect in the housing. If this happens, wait a little while until it has adjusted to the new ambient temperature before switching it on.
 - Only plug the control unit into properly installed sockets where the mains voltage corresponds to the voltage specifications on the rating plate of the control unit. Only use appropriate extension cables which comply with the technical specifications of the control unit.
 - The socket where the control unit is plugged in must be readily accessible so that it can be quickly disconnected from the mains.
 - Unwind the mains cable fully before connecting the control unit.
 - Make sure that the **control unit is always free-standing** and is not covered with a blanket or other objects.
-
1. Choose a place where you can lie comfortably and stretch out your legs, e.g. on a sofa, bed or armchair with a footrest.
 2. In order to be able to operate the control unit lying down, place it next to you on a side table or chair.
 3. Insert the mains plug of the control unit into a socket. If there isn't one nearby, use an extension cable.

B) Connecting the Cuffs to the Control Unit and Applying

If you are wearing trousers, take them off before applying the cuffs so as to ensure the treatment is effective and not unpleasant. **If you are wearing tights or compression stockings, you can keep these on.**

1. Spread out the cuffs at the point where the treatment is to be carried out.

2. The connection tubes lead to a joint connector.

IMPORTANT: The label printed on this connector must be facing upwards when it is connected so that the device works properly and the cuffs inflate upwards from the feet!

Insert the connector into the socket on the side of the control unit.

3. Open the zips of the cuffs.

4. Lie down on your back and place your legs in the cuffs.

5. Close the zips and the hook-and-loop fasteners of the cuffs. The cuffs should not be too tight around your legs (before they are inflated).

C) Switching On the Control Unit to Start Massaging

PLEASE NOTE!

- If the application becomes **unpleasant, reduce the pressure or stop** the treatment.
- If an increased pulse rate or breathing difficulties occur while using the device, stop using it. If this happens, please seek medical advice.
- **Do not walk or stand** while wearing the cuffs!
- Make sure that the **cuffs and connection tubes are not bent or twisted** so as to avoid impairing the system and damaging the air chambers!
- **Do not open the cuffs while using the device.** This could damage them.

1. Touch the  button to switch on the control unit.
2. Select mode **B (BASIC)**. In this mode, the cuffs are inflated from the feet upwards, with the pressure being maintained in the lower chambers. Any fluid in the legs is pushed towards the heart.

 In mode **A (PRO)**, the pressure "migrates" from one chamber to the next and is not maintained in the bottom chambers. Select mode **A (PRO)** only with sports applications, for example when one of the leg areas is to be left out due to an injury.

3. Make sure that all four chambers are activated and the buttons  to  on the control panel are lit. Touching the buttons will activate or deactivate the associated chamber.
If you want to leave out one area of the leg (e.g. an injured knee), deactivate the corresponding chamber on the control panel.
4. Use the  or  button to adjust the pressure in the cuffs. The pressure can be adjusted in 20 mmHg increments from 20 to 200 mmHg. However, we recommend that you do not exceed 120 mmHg.

 If a low pressure is set, the cuffs may only be fully inflated after two massage cycles.

5. Choose whether the treatment should be ended automatically after 10, 20 or 30 minutes. Touch the corresponding timer button to do this.
6. Start the treatment by pressing the  button. The control unit will then start to fill the cuffs with air and your legs will be massaged by the inflating and deflating air cushions.
7. Lie down on your back, stretch out your legs and take up a comfortable position. To ensure that the blood can circulate unhindered, it is essential to ensure that the leg in the cuff is stretched out during treatment!
8. If the pressure exerted on the legs is too powerful or unpleasant, reduce it with the  button. The ideal pressure on your legs should be powerful but not unpleasant.

 It is natural for your legs to sweat in the cuffs. This is not a cause for concern: in fact it is an intended effect which stimulates metabolism in the legs.

D) Finishing the Treatment

1. If you would like to stop the treatment early or pause it, press the  button. Otherwise the device stops automatically after the set time has elapsed.
2. Switch off the control unit entirely with the  button. NOTICE: After 10 minutes without activity, the control unit automatically switches off.
3. Wait until the cuffs have loosened. You can then open the zips and remove your legs.

E) After Treatment

- You can leave the **Sports Only from VenenWalker®** where it is for the next treatment. Be sure to remove the mains plug from the socket, however.
- If you have been sweating, clean the cuffs (see the "Cleaning" chapter).

CLEANING

PLEASE NOTE!

- Do not use any caustic or abrasive cleaning agents. These may damage the surfaces.
 - No moisture must get inside the connection tubes and the control unit.
 - Follow the care instructions on the label in the cuffs. Never wash in the washing machine!
 - Allow to dry in the air only or wipe with a cloth. Never use a hair dryer or similar device for drying! Never place the device in a tumble dryer!
-

1. Clean the cuffs with a clean, damp cloth and if necessary a mild household cleaning agent.
2. If the cuffs are used directly against the skin by more than one person, they should be cleaned using a commonly available disinfectant after use so as to avoid cross-contamination.
3. The control unit and connection tubes should only be wiped with a clean, damp cloth if necessary. You can also use a mild household cleaning agent.
4. Allow all parts to dry in the air after cleaning or use a cloth.

STORAGE

PLEASE NOTE!

- Allow all parts to dry completely before they are stored away.
-
1. Pull the connector for the connection tubes out of the control unit. Press the rest of the air out of the cuffs.
 2. Carefully fold the cuffs and connection tubes. **Avoid sharp bending, otherwise breakage points may be created on the cuffs and connection tubes after frequent use.**
 3. Store the device (in its original packaging) at temperatures of between 15 and 25°C out of the reach of children and animals. (Please refer to the "Technical Data" chapter for more details.)

TROUBLESHOOTING

If the device stops working properly, check whether you are able to correct the problem yourself based on the tips below. If this is not possible, contact the customer service department.

Do not attempt to repair an electrical device yourself!
The device may only be repaired by an authorised expert!

The control unit cannot be switched on.

The control unit is not connected to the mains.

- ▶ Insert the mains plug into the socket.

There is no mains voltage.

- ▶ Try another mains socket.

The mains cable is damaged.

- ▶ Check the mains cable. If it is damaged, contact the customer service department.

Only one cuff or no cuff at all is inflated or the pressure is very low even when it is set to be high.

The connector for the connection tubes is not connected correctly.

- ▶ Slide the connector into the socket on the side of the control unit as far as it will go.

Air escapes between the connection tubes and cuffs.

- ▶ Check on the cuffs whether the connection tubes are plugged correctly on the tube connectors. If they are not, simultaneously twist and slide the end of the tube onto the tube connector. Repeat this for the other connection tubes if necessary.

The connection tubes are not tight.

- ▶ Check whether the connection tubes are porous, have any holes or other damage and contact the customer service department if necessary.

The air input is obstructed.

- ▶ Spread out the cuffs and connection tubes. Make sure that the connection tubes and cuffs are not bent, twisted or compressed.

The cuff is not tight.

- ▶ Check the cuff for any tears, holes or other damage and contact the customer service department if necessary.

The cuffs do not deflate.

The air flow is obstructed.

- ▶ Make sure that the connection tubes are not bent, twisted or compressed. Pull the connector for the connection tubes out of the control unit and press the rest of the air out of the cuffs.

The control unit is unusually loud.

The control unit is placed on an unsuitable surface.

- ▶ Place the control unit on a level and firm surface.

TECHNICAL DATA

Product name:	Sports Only from VenenWalker®
Model number:	L-1011
Reference article no. for customer service:	14571
Mains voltage:	230 V ~
Mains frequency:	50 Hz
Current consumption:	0.2 A
Power:	50 W
Cuff pressure:	max. 200 mmHg
Weight of control unit:	2600 g
Fuse:	T2 A H250 V (5 x 20 mm)
Electric shock protection class:	Type BF 
Fluid ingress protection class:	IP21
Protection class:	II
Temperature range:	Operation: +5 – +35°C Transport: -20 – +55°C
Humidity:	Operation: 10 – 90% Storage and transport: max. 90%
Air pressure:	Operation: up to altitude of 2000 m, 80 – 106 kPa Storage and transport: 80 – 106 kPa
Air volume:	max. 16 – 18 l/min
Dimensions:	- Control unit: 285 mm (L) x 175 mm (W) x 130 mm (H) - Cuffs: 85 cm (L), thigh circumference approx. 70 cm, calf circumference approx. 42 cm - Connection tubes: 160 cm
Cuff material:	Nylon (outside and inside)
Pump housing material:	ABS/PA-765; UL classification of flammability: V-0
Service life of the system:	5 years
Serial number:	GEGM14240223 (GEGM14 = internal number, 24 = 2024, 02 = February, 23 = day)
ID of operating instructions:	Z 14571 M GM V1 0524 md/as
Date of operating instructions:	03.05.2024

CUSTOMER SERVICE

If you have any questions about **Sports Only from VenenWalker®**, contact us via our website **www.venenwalker.com**

or by phone on +49 (0) 38851 314337* (Mon – Fri, 8 am – midday)***

Contact

E-mail: beratung@venenwalker.com

MORE INFORMATION

Sports Only from VenenWalker® is manufactured by:

XIAMEN SENYANG CO.,LTD

Address: No 189-1, Tongji South Road, Jimei District, Xiamen, 361000, Fujian, P.R. China

E-mail: vip@senyangcare.com

Web: www.medicalxm.com

EU Representative:

SUNGO Europe B.V.

Address: Fascinato Boulevard 522, Unit 1.7, 2909VA Capelle aan den IJssel, Netherlands

Tel: +31(0)10 3034500; +31(0)2021 11106

Importer:

Globalmind Consumer Electronics GmbH

Address: An der Strusbek 50, 22926 Ahrensburg, Germany

Tel: +49 40 7392140

Fax: +49 40 73921418

E-mail: beratung@venenwalker.com

Distributor:

DS Produkte GmbH, Stormarnring 14, 22145 Stapelfeld, Germany

*** Calls to German landlines are subject to your provider's charges.

DISPOSAL

Dispose of the packaging material in an environmentally friendly manner so that it can be recycled.

The adjacent symbol (crossed-out bin underlined) means that used devices should not be disposed of in the domestic waste, but should instead be taken to special collection and return systems.

Owners of used devices from private households can take them free of charge to the collection points of the public waste disposal authorities or the disposal points which are set up by manufacturers and distributors under the Electrical Equipment Act so that they can be disposed of in an environmentally friendly way and valuable raw materials can be recovered. If they are not disposed of properly, toxic ingredients may escape into the environment and cause harm to the health of humans, animals and plants. Businesses which market electrical and electronic devices are also obliged to take back used devices.

Waste prevention makes an even more valuable contribution to protecting the environment. This means if possible, as well as continuing to use or repair a device yourself, passing it on to another user also represents an ecologically beneficial alternative to disposal.

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EMC Manufacturer's Declaration

Model PT1003 electromagnetic emissions—manufacturer's declaration		
The model PT1003 is intended for use in the electromagnetic environment specified below. The customer or the user of the model PT1003 should assure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment—guidance
RF emissions CISPR 11	Group 1	The model PT1003 is suitable for use in all establishments, including domestic establishments and those directly connected to the public low voltage power supply network that supplies building used for domestic purposes.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations / flicker emissions IEC 61000-3-3	Not applicable	

Guidance and manufacturer's declaration—electromagnetic immunity			
The model PT1003 is intended for use in the electromagnetic environment specified below. The customer or the user of the model PT1003 should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment—guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 2, 4, and 6 kV contact ± 2, 4 and 8 kV air	± 2, 4, and 6 kV contact ± 2, 4 and 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrostatic fast transient / burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input / output lines	± 2 kV Not applicable	Mains power quality should be that of a typical home use location.
Surge IEC 61000-4-5	± 0.5 and 1 kV line(s) to line(s) ± 2 kV line(s) to earth	± 0.5 and 1 kV differential mode Not applicable, no ground wire	Mains power quality should be that of a typical home use location.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5 % U_T (>95% dip in U_T) for 0.5 cycles 40 % U_T (60% dip in U_T) for 6 cycles 70 % U_T (30% dip in U_T) for 30 cycles <5 % U_T (>95% dip in U_T) for 5 s	<5 % U_T (>95% dip in U_T) for 0.5 cycles 40 % U_T (60% dip in U_T) for 6 cycles 70 % U_T (30% dip in U_T) for 30 cycles <5 % U_T (>95% dip in U_T) for 5 s	Mains power quality should be that of a typical home use location. If the user of the model IC-1545-DL requires continued operation during mains power interruptions, it is recommended that the model IC-1545-DL be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	The power frequency magnetic fields should be at the levels found in a typical home use location.
NOTE: U_T is the AC mains voltage prior to application of the test level.			

Guidance and manufacturer's declaration — electromagnetic immunity

The model PT1003 is intended for use in the electromagnetic environment specified below. The customer or the user of the model PT1003 should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment—guidance
---------------	----------------------	------------------	--------------------------------------

Conducted RF IEC 61000-4-6	3 Vrms 15 kHz to 80 MHz	3 Vrms 3 V/m	Portable and mobile RF communications equipment should be used no closer to any part of the model PT1003, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance 80 MHz to 800 MHz 800 MHz to 2.5 GHz where <i>P</i> is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and <i>d</i> is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,» should be less than the compliance level in each frequency range.» Interference may occur in the vicinity of equipment marked with the following symbol: 
Note 1: At 80 MHz and 800 MHz, the higher frequency range applies. Note 2: These guidelines may not apply to all situations. Electromagnetic propagations is affected by absorption and reflections from structures, objects and people.			
»Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. if the measured field strength in the location in which the model PT1003 is used exceeds the applicable RF compliance level above, the model PT1003 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the model PT1003. »Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.			

Recommended separation distances between portable and mobile RF communications equipment and the Model PT1003

The model PT1003 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer of user of the model PT1003 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the PT1003 as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	$d = 1.2\sqrt{P}$	$d = 1.2\sqrt{P}$	$d = 1.2\sqrt{P}$
	150 kHz to 80 MHz	80 MHz to 800 MHz	800 MHz to 2.5 GHz
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance *d* in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where *p* is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note 1: At 80 MHz, the separation distance for the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

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All rights reserved.
DS Produkte GmbH, Stormarnring 14, 22145 Stapelfeld, Germania