VITALmaxx

Shiatsu-Massagegerät

15483

DE

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EN

Operating instructions
from page 11

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LIST OF CONTENTS

Symbols	11
Signal Words	11
Further Explanations	11
Items Supplied	13
At A Glance	13
Intended Use	13
When May the Device Not Be Used?	14
Safety Notices	14
Charging Up the Battery	16
Use	16
Cleaning and Storage	17
Technical Data	17
Disposal	18

IMPORTANT NOTICES! PLEASE KEEP FOR REFERENCE!

SYMBOLS

Danger symbols: These symbols indicate possible dangers. Read the associated safety notices carefully and follow them.

Supplementary information

Read the operating instructions before use!

Direct current

Alternating current

SIGNAL WORDS

Signal words indicate hazards if the associated notices are not followed.




WARNING – medium risk, may have serious injury or death

NOTICE – may result in risk of damage to material

FURTHER EXPLANATIONS

Key numbers are shown as follows: **(1)**

Picture references are shown as follows: **(Picture A)**

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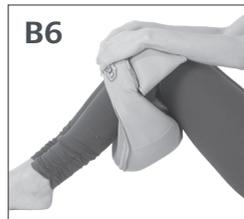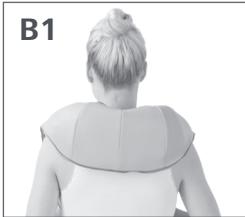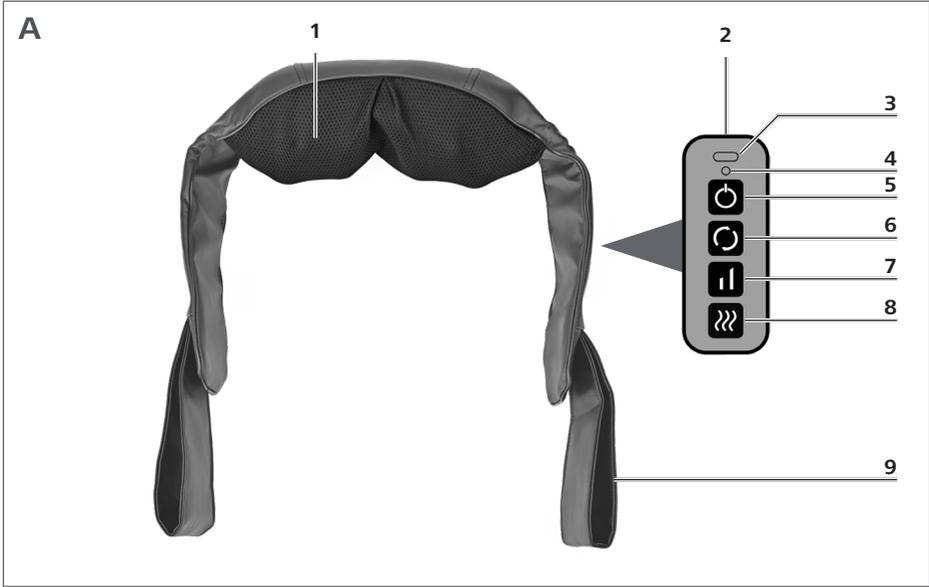Information About the Operating Instructions

Before using the device for the first time, please read through these operating instructions carefully and keep them for future reference and other users. They form an integral part of the device. The manufacturer and importer do not accept any liability if the information in these operating instructions is not complied with.

If you have any questions about the device and about spare parts/accessories, please contact the customer service department via our website: www.ds-group.de/kundenservice

ITEMS SUPPLIED

- 1 x massage device with removable cover
- 1 x Charging cable (USB-A to USB-C)
- 1 x operating instructions

Check the scope of delivery for completeness and the components for transport damage. If you find any damage, do not use the product but contact our customer service department.

Remove any possible films, stickers or transport protection from the product. **Never remove the rating plate and any possible warnings!**

A USB mains adapter is required for charging at the plug socket. This is not supplied with the device.

AT A GLANCE

(Picture A)

- 1 Massage balls, heating and red light element
- 2 Control panel
- 3 Charging socket
- 4 Control lamp:
 - red (flashing): battery is charging / battery needs to be charged
 - green (constant): battery is charged / device is switched on (**without** heat)
 - red (constant): device is switched on (**with** heat)

- 5 button: Switch device on/off
- 6 button: Change direction of the massage balls manually
- 7 button: Change speed (2 levels)
- 8 button: Switch heat and red light on/off
- 9 Hand loops

INTENDED USE

- This device is intended to be used for loosening the neck and back muscles and delivers relaxation. The device may also be used to massage abdominal, buttock and leg muscles. It is not a replacement for medically required treatment.
- The device is for personal use only and is not intended for commercial applications.
- The device should only be used as described in these instructions. Any other use is deemed to be improper.
- The following are excluded from the warranty: all defects caused by improper handling, damage or unauthorised attempts at repair. This also applies to normal wear and tear.

WHEN MAY THE DEVICE NOT BE USED?

- The following people may use the device **only following consultation with a doctor**:
 - people who have no feeling in their neck, back, buttocks and / or thighs / lower legs
 - people who have back complaints or spinal disorders
 - people with circulation problems
 - women who are pregnant or have just given birth
 - people with a pacemaker
- people who have undergone an operation for which the treatment is not yet complete.
- **Only use the device if you have no physical illnesses or ailments!** If you are unsure whether the device may be used, consult a doctor.
- Do not use the device near to or on injuries, swelling, skin irritations, infected skin, warts, fungal infections, abnormal changes in the skin or body parts with no feeling.

SAFETY NOTICES

WARNING: Note all safety notices, instructions, illustrations and technical data provided with this device.

Failure to follow the safety notices and instructions may result in electric shock, fire and / or injuries.

- This device may be used by **children** aged 8 and over and by people with reduced physical, sensory or mental abilities or a lack of experience and/or knowledge if they are supervised or have received instruction on how to use the device safely and have understood the dangers resulting from failure to comply with the relevant safety precautions.
- **Children** must **not** play with the device.
- Cleaning and **user maintenance** must not be carried out by **children**, unless they are older than 8 years of age and are supervised.
- When the heat function is switched on, the device has a hot surface. People who due to illness or for other reasons are not sensitive to heat must be careful when using the heat function.
- The device must always be disconnected from the mains power when it is left unattended and before cleaning it.

- Do not make any modifications to the device. If the power cord of this appliance is damaged, it must be replaced by a special power cord available from the manufacturer or its service agent.

WARNING – Danger of Electric Shock

- Never immerse the device and the charging cable in water or other liquids and ensure that these components cannot fall into water or become wet.
- If the device falls into water while it is connected to the mains power, pull the USB mains adapter out of the plug socket immediately! **Only after this** should you pull the device out of the water. Have it checked in a specialist workshop before using it again.
- Never touch the device and the charging cable with wet hands when these components are connected to the mains power.

WARNING – Danger of Fire

- Do not cover the device during operation or charging in order to prevent it from catching fire.
- Do not use the device under a blanket or similar item or in the vicinity of combustible material.

WARNING – Danger from Rechargeable Batteries

- The device has a built-in rechargeable battery. It must not be thrown into a fire or exposed to high temperatures of over +40°C. Do not expose it to direct sunlight for a prolonged period of time. Do not disassemble the battery! Do not place any heavy objects on the device. The battery could catch fire or explode. Do not store the device with a built-in battery in a vehicle where it might be exposed to extreme temperatures.

- If the battery leaks, do not allow your skin, eyes and mucous membranes to come into contact with the battery acid. Wear protective gloves if necessary. If you touch battery acid, rinse the affected area with abundant fresh water immediately and contact a doctor as soon as possible. There is a danger of burns from battery Acid!
- Do not use the device if the battery is leaking.
- Protect the device from persistent moisture. Do not use, charge or store it in rooms with high humidity.

NOTICE – Risk of Damage to Material and Property

- Only connect the mains adapter to a plug socket that matches the technical data of the device.
- Only use suitable extension cables whose technical data is the same as that of the device.
- In the event of damage / faults (e.g. unusual noises, odour or smoke development), switch the device off immediately and disconnect it from the power supply.

CHARGING UP THE BATTERY

Please Note!

- Do not charge on a computer because the high charging current could damage the computer.
- **It is essential that you charge up the battery fully before you first use the device in order to achieve the maximum charging capacity.** If the battery is fully discharged during use, it must be fully charged up again before it is next used.

If the device cannot be switched on or the control lamp (4) flashes red, the battery needs to be charged:

1. Plug the charging cable into the charging socket (3) and then connect it to a power source with a USB interface (preferably via a USB mains adapter to the plug socket). The control lamp ...
 - flashes red: battery is charging
 - lights up green: battery is fully charged
2. Disconnect the charging cable from the power source and then from the charging socket.

USE

Please Note!

- For reasons of hygiene, avoid direct contact of the skin with the device wherever possible. It should not come into contact with skin that has had cream or oil applied to it.
- Do not insert any fingers or other body parts into the device. Do not open the zip during use.
- No pressure should be exerted on the device as this will slow down or stop the massage and the device could be damaged. Do not stand or lie on the device.

- Switch off the device immediately if you experience pain or feel unwell during use.
- A tingling sensation/itching during or after the massage is normal due to the circulation-boosting effect.
- Do not use the device for longer than 30 minutes without stopping.
- Do not massage one group of muscles for longer than 15 minutes.
- Do **not** fall asleep while you are using the device!

Do not use the device before going to bed. The massage has a stimulating effect, which may lead to problems getting to sleep. We recommend that you wait an hour after mealtimes before using the device.

Areas of Application

The device is designed in particular for use on the upper back area (in particular the neck; **Picture B1 + B2**). But it can also be used to massage abdominal, buttock and leg muscles (**Picture B3–B6**).

Operating the Device

1. Make sure the device is disconnected from the power supply.
2. Place the massage surface (1) of the device onto the part of the body which is to be massaged. For a shoulder or back massage, place the straps of the device around your shoulders and fix them in place by putting your arms through the loops (9) (**Picture B7**) and thus holding the device in the correct position.
3. **Switch on:** Hold down the  button (5) for approx. 2 seconds. The red light and heat function is automatically switched on at the same time. The massage balls massage either clockwise or anticlockwise. After approx. 1 minute the direction changes.

- **Change direction manually:** Briefly press the  button (6).
 - **To change speed (2 levels):** Briefly press the  button (7).
 - **Switch heat on/off:** Briefly press the  button (8).
4. After 15 minutes the device automatically switches off. To switch it off sooner, hold down the  button for approx. 2 seconds.
 5. Allow the device to cool down if the heat function was activated before it is stored away or cleaned (see the “Cleaning and Storage” chapter).

CLEANING AND STORAGE

Please Note!

- Do not immerse the device or the cable in water!
 - The device may not be connected to the mains for cleaning!
 - Do not use any caustic or abrasive cleaning agents for cleaning. These can damage the surface.
-
- Wipe down the device just with a slightly damp cloth if necessary. Then wipe it down fully with a dry cloth.
 - Keep the device in its respective packaging in a dry place which is protected from sunlight and is not accessible to children and animals.
 - If unused for a prolonged period of time, recharge the battery completely every three months.

TECHNICAL DATA

Article number:	15483
ID operating instructions:	Z 15483 M DS V1.1 0325 rr
Model number:	C37-MP085-R
Voltage supply:	7.4V DC (Li-Ion battery 2000 mAh)
Power:	max. 24W
Battery charge time:	approx. 3–4 hours
Battery operating time:	max. 45 minutes (with fully charged battery, depending on mode)
Timer:	15 minutes
Mains adapter required:	5V and 2A

DISPOSAL

Dispose of the packaging material in an environmentally friendly manner so that it can be recycled.

The adjacent symbol (crossed-out bin underlined) means that used devices should not be disposed of in the domestic waste, but should instead be taken to special collection and return systems.

Owners of used devices must remove from the used device any used batteries or rechargeable batteries which are not encapsulated by the used device and can be removed in a non-destructive way and dispose of them separately (see also the section on battery disposal).

Owners of used devices from private households can take them free of charge to the collection points of the public waste disposal authorities or the disposal points which are set up by manufacturers and distributors under the Electrical Equipment Act so that they can be disposed of in an environmentally friendly way and valuable raw materials can be recovered. If they are not disposed of properly, toxic ingredients may escape into the environment and cause harm to the health of humans, animals and plants. Businesses which market electrical and electronic devices are also obliged to take back used devices.

(Rechargeable) batteries must not be disposed of along with the household waste. Consumers are legally obliged to dispose of (rechargeable) batteries separately.

(Rechargeable) batteries can be handed in free of charge at a collection point of the local authority/district or in retailers so that they can be disposed of in an environmentally friendly way and valuable raw materials can be recovered. If they are not disposed of properly, toxic ingredients may escape into the environment and cause harm to the health of humans, animals and plants.

Only dispose of (rechargeable) batteries when they are discharged. If possible, use rechargeable batteries rather than single-use batteries.

In the case of (rechargeable) batteries containing lithium, tape over the poles before you dispose of them to prevent a short circuit. A short circuit may cause a fire or explosion.

Waste prevention makes an even more valuable contribution to protecting the environment. This means if possible, as well as continuing to use or repair a device yourself, passing it on to another user also represents an ecologically beneficial alternative to disposal.

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☎ +49 38851 314650*

kundenservice@dspro.de

* Calls to German landlines are subject to
your provider's charges.

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All rights reserved.
DS Produkte GmbH